Jordan CLEAN PREMIUM

Electrical Toothbrush

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5 V DC-lader, 800 mA

Valmistaja:

WILFA AS
Industriveien 25,
1481 Hagan
Norway

JORDAN CLEAN ELECTRIC TOOTHBRUSH

The Jordan Clean series is a new generation of electric toothbrushes that through SoftClean™ Technology delivers the right combination of oral hygiene features. Inspired by a Scandinavian design approach for simplicity and user-friendliness, the series gives an effective clean from an electric toothbrush and the ease of use of a manual one.

SAFETY

Please note the following instructions when using the appliance:

- This appliance is intended for domestic use only.
- This appliance is designed solely for cleaning and caring for human teeth. All other uses are expressly prohibited.
- Do not use if the appliance or attachments is visibly damaged.
- When the appliance is being charged, on no account should this or the charging unit come into contact with water or any other liquid.
- Keep the appliance out of direct sunlight, excessive heat and open flames.
- Do not drop the appliance.
- Keep the appliance out of the reach of children.
- Do not open the appliance under any circumstances.
- No warranty claims are accepted for damage caused by improper handling.
- Cleaning and user maintenance shall not be made by children without supervision.
- This appliance can be used by children aged from 8 years and above. And persons with reduced physical, sensory or mental capabilities, or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.

SPECIAL FEATURES:

- **Softclean™ technology** gives you a unique combination of highly effective bristles and a smaller head for effective cleaning in a gentler way. Each head is divided into three cleaning zones with specialized efficiencies per zone.
- **Built-in high capacity Lithium battery** for long use. It can run 120min which gives 1 month based on 2 times brushing a day and 2 minutes brushing each time.
- **Quick charge** with at least 20 minutes you can use it for nearly a week.
- **Recharging the brush** takes 4-6 hours for a full charge depending on how much left on battery. Automatic cut off when full charged.
- **Detachable USB cable** gives a hassle free charging environment. Toothbrush can be charged into any USB socket.
- **3 modes** with specially designed Whitening

mode with oscillating speed at 10 000 rpm, Sensitive mode with oscillating speed 8000 rpm and Clean mode at 9000 rpm. Hold the switch for 1s for selecting the desired Speed.

- **Pressure sensor** when excessive pressure is applied when brushing, the pressure sensor indicator will appear on LED display and the oscillating speed will be reduced automatically to remind you to lower your pressure and also protecting your teeth and gum.
- **Travel case** is included for convenient use with space for 2 spare toothbrush heads.
- **2-min timer** with click signal for each 30 sec brushing time for reminding you to move to different teeth area. It will click 2 times at 30/60/90 sec. and click 4 times at 2 minutes (recommended brushing time). The brush will stop automatically after 10 minutes for protecting your teeth and preserve battery capacity by accidental switch-on when travelling.
- **Speed setting memory.** Always memorize the last Speed for next use.
- **Low power indication** when battery capacity goes down to about 20%, the indicator blinks quickly to remind you to recharge the brush in order to keep full power.
- **Over-charge and Over-discharge Protection** circuit is equipped to preserve the battery life time. When full charge the toothbrush will automatic cut off. When battery is too low the toothbrush will automatic stop.

SPECIFICATIONS

Adaptor:100-240V~ 50/60Hz 4W. Rechargeable power unit is 3.7Volts. Maximum operation time is approximately 120 minutes/1 month when fully charged. When rechargeable batteries are running low, the indicator will blink to tell you to recharge the unit immediately. If you continue to use the unit for some time without recharge, the unit may stop automatically. This is our battery protection feature to preserve the batteries life time. You need to recharge immediately.

CONNECTING

Insert the USB cable plug into the adaptor. Then the adaptor can be connected to any voltage outlet - even to shaver sockets - between 100-240 volts AC (~). Connect the other end of USB cable to the charging base. The indicator will blink with light/shade(breathing mode) alternately, when the power unit is inserted on the plugged-in charging base. It will light up WHITE when batteries are fully charged (4~6 hours). Besides, you can also connect the USB cable to the computer directly for charging. The power unit may warm up somewhat during recharging. This has no significance on the unit's operation. Over-charging is impossible.

BEFORE USE

Before using your Jordan Clean toothbrush for the first time, charge the unit full. After the power unit has been fully charged, attach the brush head to the power unit. The brush head must face the same direction as the switch after twist in. Jordan Clean has only one switch by pressing for the desired functions. In case you put the power unit into the charging base when it is working, it will automatically stop.

DISPOSAL OF RECHARGEABLE BATTERY

This toothbrush is using Lithium rechargeable battery, after the end of the toothbrush life time, you need to deliver toothbrush for proper disposal according to your local authority instruction for environmental protection disposal. The battery must be removed from the appliance before it is discarded. The appliance must be disconnected from the mainspower

supply before removing the battery. The battery must be disposed of safely. Use a flathead screwdriver to remove the internal body, then remove the battery and dispose of it via a return and collection system or contact the retailer from whom the product was purchased. They operate a take-back system for safe recycling.

PROTECTING THE ENVIRONMENT

Correct Disposal of this product.

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

POWER SUPPLY

Adaptor

Input:
100-240 V~, 50/60 Hz, 0,2 A


Output:
5 V 800 mA

Appliance

Li-polymer battery 3.7V, charged by DC5V, 800mA

Produced by:

WILFA AS
Industriveien 25,
1481 Hagan
Norway

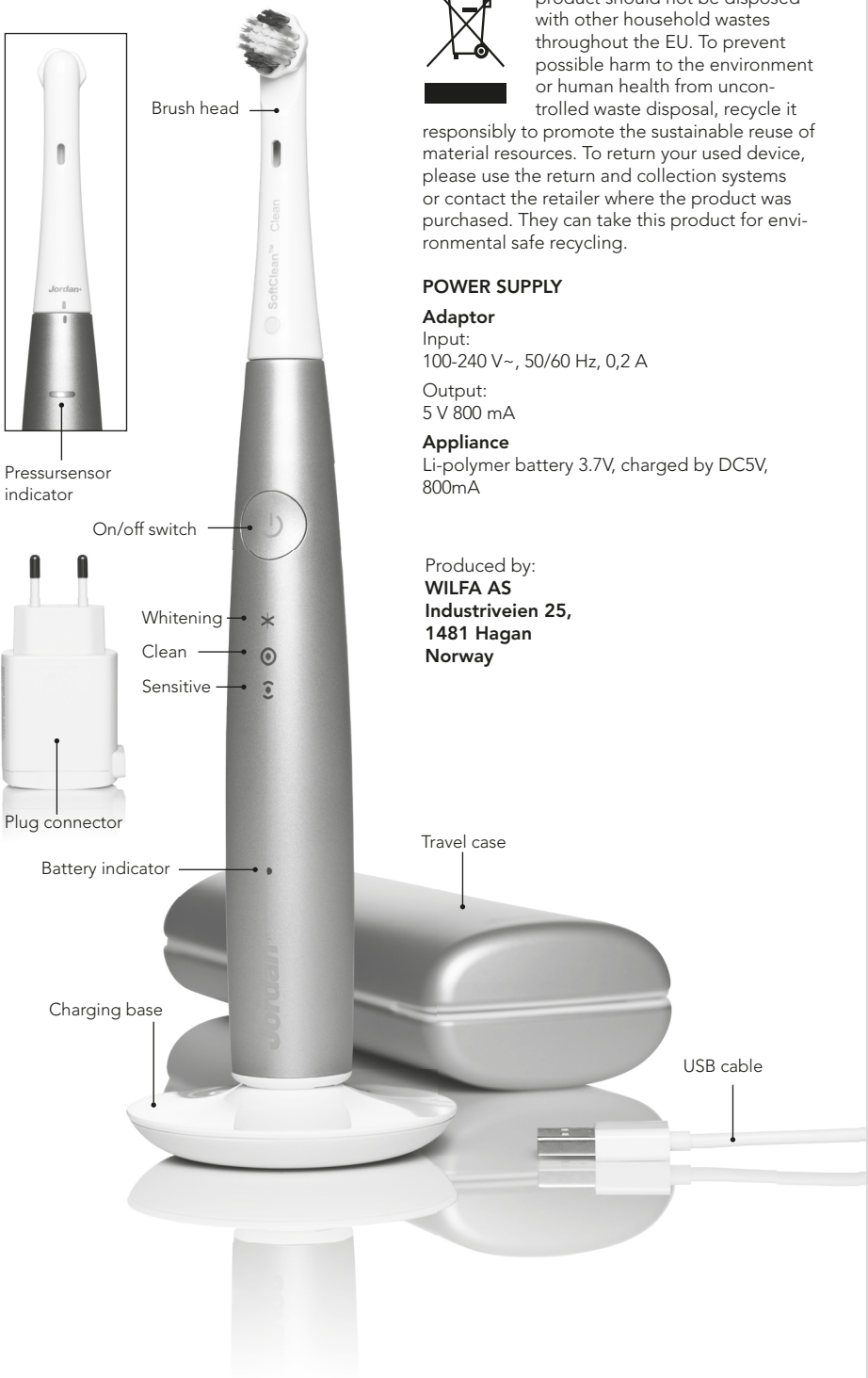Markedsføres i Norge av:
WILFA AS
Industriveien 25
Postboks146, Hagan
1481 Hagan, Norway
Tlf: 67 06 33 00
support@wilfa.com
www.wilfa.com

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WILFA AB
Oy Wilfa Suomi Ab